

Cobolt FIC-07

Mounting plate for laser-to-fiber coupling - no fan

Cobolt model number	FIC-07
Thermal resistance	< 0.25 K/W
Material	Aluminium
FIC-07 is suitable for use with all Cobolt lasers in the following models:	Cobolt 04-01, 06-01, 08-01 Series
Dimensions [mm] [inches]	216 x 118 8.5 x 4.6
Weight	1,8 kg



The option FIC-07 is recommended for Cobolt 04-01, 06-01, 08-01 Series. The option consists of a non-fanned heatsink and a coupling block to mount the fiber coupler onto.

Fiber coupling options for either multimode fiber (MM) or single-mode polarization maintaining (SM/PM) fiber can be ordered separately. The FIC-07 can be bolted to a standard metric or imperial optical table.



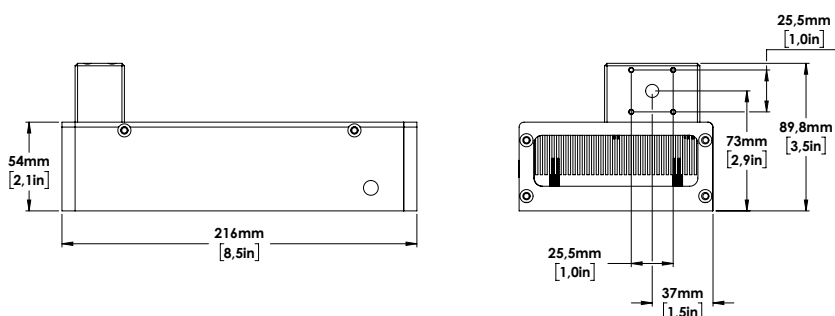


Illustration of possible configuration. Laser and fiber coupling option should be ordered separately.

HÜBNER Photonics | Coherence Matters.

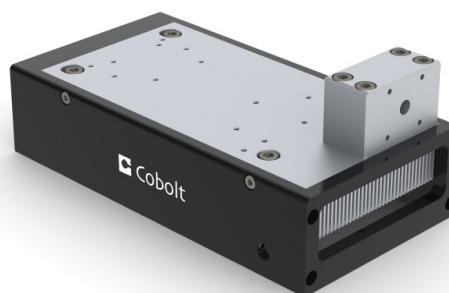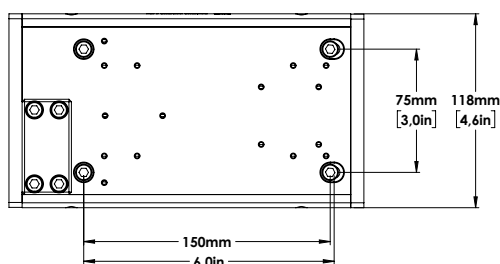